

# Ka la a'es'ao'āk tapen Rapid Antigen Test

Fas'āk tē ne la nā hūn se 'on la a'sok tapen ne test te'is la tū'āk se' 'e kəinag brand tūtū 'atakoa. Rapid antigen testing kits ma'oi tāla ʔaupir se fas'āk tē ne fā' sin 'e lopo hūn se 'on la a'sok ne tēkāe te'is, ka ne 'eagke'ea 'e av 'atakoa. Figalelei ma ʔaupir se fask'āk tē ne kampanē ne rēam kit heta (manufacturer) fā' sin.

## Teak ne ta sample het




Hot'ākim ta nasal swab het 'e laloag tag het ta (pouch).



Tem'āk swab heta se laloag ne is hap 'eseat la hāiatəuəg ma sentimīta 2-3 'e utut te palag ne 'ou isu.



Takmūl'āk a'liəm (5 times) swab heta 'e laloag ne 'ou is hapat 'Āe la a'es'ao'ākia swab 'eseat ma ne kotā a'es'ao'āk sio la a'sok hoi'ākia 'amnāk te'is se 'ou is hapat te tore.

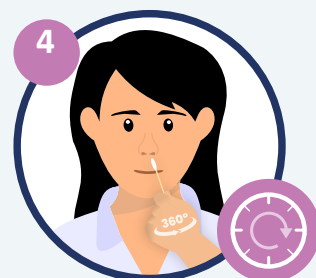

Səkiar se fas'āk tē ne leum 'e kes he ta (kit box instructions) la a'fūmou'ākia 'on roa ne avat ne 'āe la ʔari ma kotā 'āe pō la hatia roggro ne mou se səkiar tē ta (test result) 'On roa ne av ta kop ma la tū'āk sio ka ne tae ke 'e kit het ne 'āe pō.

## Səkiar ne sample



Sis ma 'ut'āk aluminium foil seal ne nā sin se rere ne extraction vial, ne hāe sin extraction buffer.




Nā swab heta se laloag ne extraction vial. Takmūl'āk a'lelei swab heta, mea'mea' se' ma ka la takmūl'āk a'liəm.




Takmūl'āk swab heta ma has'āk se vakvak ne jou heta (vial) la 'ut'ākia swab heta, has'āk se vakvak heta la hō ʔənut ne tae 'e swab heta la vil 'atakoa se laloag jou heta (vial).



'I'om vial heta 'e popou het ne nā sin ma has la mou se rer ne vial heta. Suə a'lelei ʔau ma rūrū fān ne jou heta (tube).

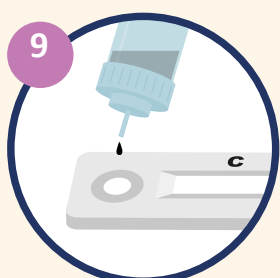

Hā'fūreū'āk jou (vial) heta se avat ne sample heta tae 'e rer ne sample well heta. Hō a'riri' jou (vial) heta. A'sok la ʔən ne sample ta la jer a'riri' se laloag ne sample well heta la hele' se avat ne jeru la 'atakoa'āk fakma se fas'āk tē ne fā' sio 'e kes (kit) heta.

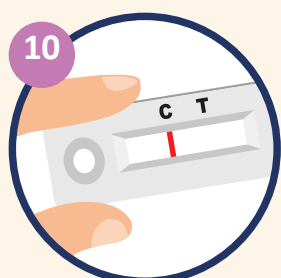

Figalelei ma səkiar a'lelei se fas'āk tē ne fā' 'e kes (kit) heta ma ʔaupir se av het ne fə'i ma kotā 'āe 'io se roggro ne mou se 'ou səkiar tē (test) ta.

## Hat ne 'ou roggro ne 'ou səkiar tē (test) ta


### Negative

Sah 'esea het 'e for ne maf ne puk te'is C kel'āk ne 'āe negative ma kat pō 'e ra 'e 'af'af ta.




### Positive

Sah he ruə, sahat 'e for ne maf ne puk te'is C ma sahat 'e for ne maf ne puk te'is T, kel'āk ne 'āe positive ma pō 'e 'af'af ta.

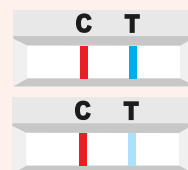

### Roggo Həisiriəg - Invalid Result

Kat 'es sah het ra no' 'e for ne maf ne puk te'is C kel'āk ne test ta həisiriəg ne invalid.

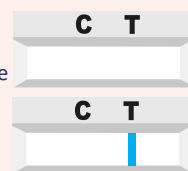

## Rak'āk ne 'ou roggro

'Āe kop la rak'āk roggro ne mou se 'ou test ta. Tē te'is la pō se' ma la a'sok 'e My Covid Record.

Log in se [mycovidrecord.health.nz](https://mycovidrecord.health.nz) ne he' se **0800 222 478** ma has se vəgat te 'on **3**.

Nōnō ka 'āe tae 'e Close Contact Exemption Scheme, 'āe noj la rak'āk se 'ou ut garueag ta 'ou roggro ne test ta.

